



cql-7---edico

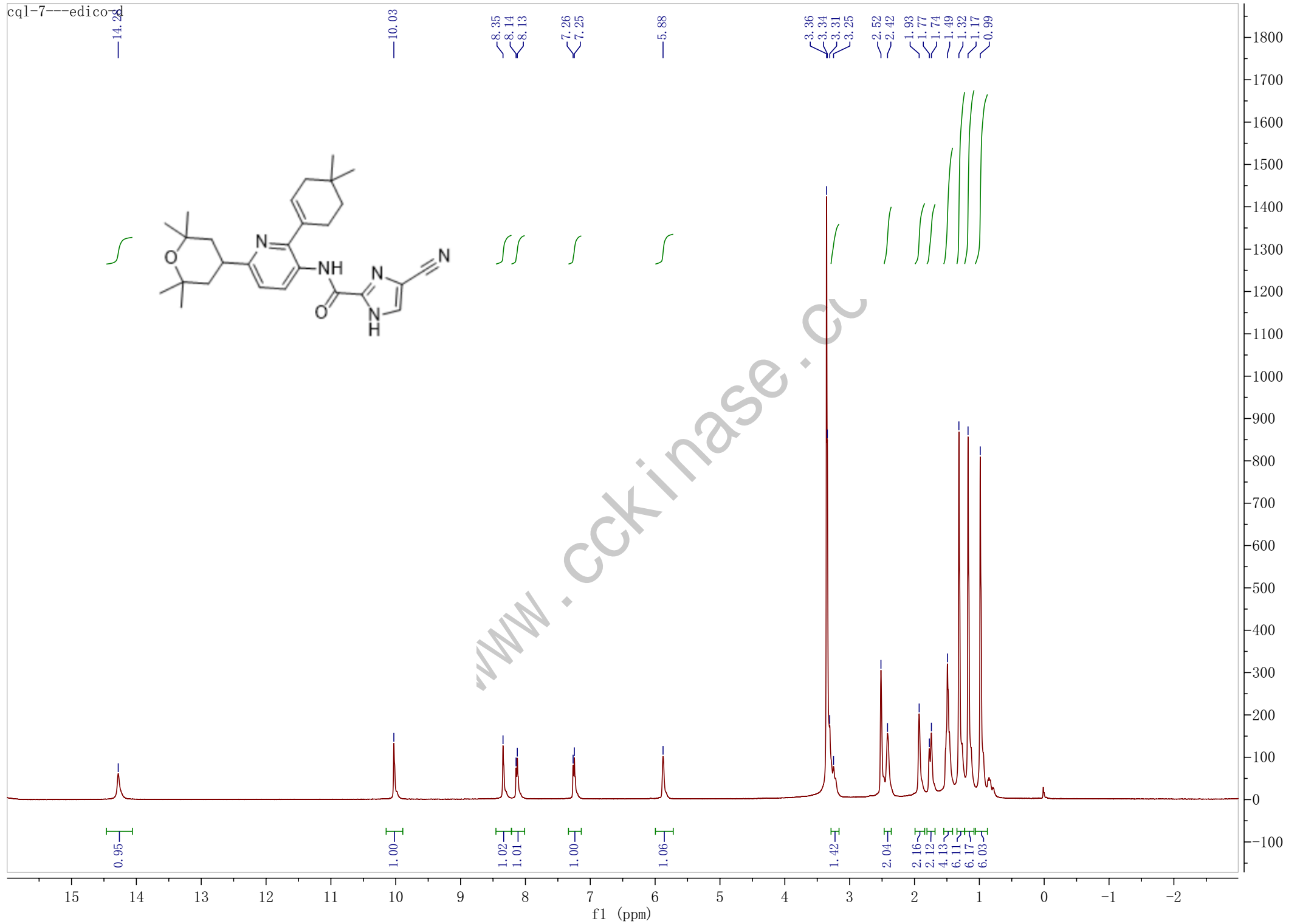

==== Shimadzu LabSolutions Analysis Report ====

Sample Name : EDICOTINIB
 Vial# : 1-10
 Date Acquired : 2020/7/3 14:59:49
 Date Processed : 2020/7/3 15:23:58
 Data Filename : D:\DATA\20200615\EDICOTINIB.lcd
 Method Filename : C:\LabSolutions\Data\Project1\LCMS-WT1.lcm
 Batch Filename : D:\DATA\20200615\20200615.lcb
 FLOW RATE: 1ML/MIN
 COLUMN: AGILENT ECLIPSE XDB-C18 150*4.6MM 5UM
 COLUME TEMPERATURE: 25
 GRADIENT PROGRAMME:
 TIME(MIN) MOBILE PHASE-A(H2O+0.1%HCOOH) MOBILE PHASE-B(CH3CN+0.1%HCOOH)
 0 90 10
 5 0 100
 20 0 100

mAU

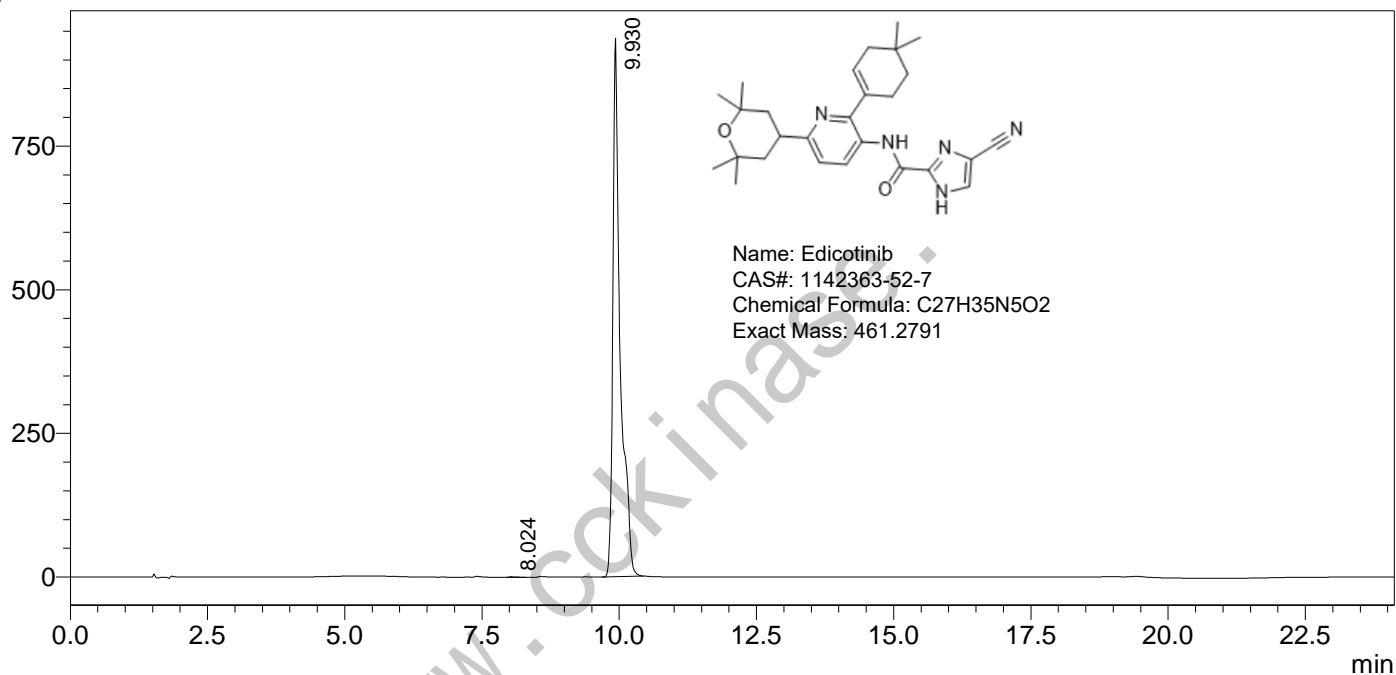

PDA Ch1 280nm

Peak#	Ret. Time	Height	Height%	Area	Area%	Name
1	8.024	948	0.101	6382	0.073	
2	9.930	937680	99.899	8773066	99.927	
总计		938628	100.000	8779448	100.000	

==== Shimadzu LabSolutions 分析报告 ====

样品名 : EDICOTINIB
分析者 : System Administrator
分析日期/时间 : 2020/7/3 14:59:49
样品类型 : 未知
样品瓶架号 : 1
样品瓶号 : 10
进样体积 : 1
数据文件 : EDICOTINIB.lcd
方法文件 : LCMS-WT1.lcm
处理者 : System Administrator
处理日期/时间 : 2020/7/3 15:23:58

#: 1 保留时间: 9.967(扫描#: 599)

峰数: 977

质谱: 单图 9.967(599) 基峰: 462.15(8390521)

背景: 无 段 1 - 事件 1

